

Product Change Notification


Product Change #: **PCN326**

Posting Date: Jan 27th

Description of Change: Nordson Medical part number FCA026-001 stopcock tubing assembly is going through a change with respect to the stopcock housing used in the assembly. These housings will now be annealed to reduce the high molded-in stresses and thus prevent environmental stress cracking of these housings that can occur over a period of time.

It is important to note that the raw materials and finished good part numbers will not change.

For your reference, the specific notification number of this change is PCN326.

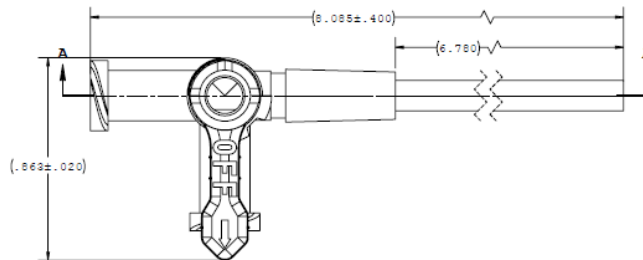

Reason: Nordson Medical part number FCA026-001 stopcock tubing assembly can crack at the tubing bond location. These cracks are determined to be due to high molded-in stresses within the polycarbonate housing, eventually relieved when a solvent is used to bond the tubing. Root cause and problem solving analyses has proven that annealing the stopcock housings is the solution.

Parts Affected: FCA026-001

Effective Date: Jan 15th, 2015

Implementation Procedure:
All current finished goods need to be scrapped. This will be a rolling change.