



SPAULDING SERIES™ **Aseptic Disconnects**

Providing a leak-free connection, the Spaulding Series™ Aseptic Disconnects eliminate the need for pinch clamps and tube welders. Smooth disconnects are achievable with the easy-to-use design, and valuable fluids are protected by valves that automatically shut off flow when the fittings are disconnected. Available in USP Class VI, ADCF Polysulfone, validation and test results confirm their suitability in pre-sterilized systems and are also gamma radiation stable. These connectors are interchangeable with similar plastic connectors.


FEATURES & BENEFITS

- Valves that are automatically shut off flow at disconnection
- Interchangeable with other plastic connectors
- Plugs and caps available to maintain fluid path integrity
- Made with BSE/TSE-free materials
- Validation and test results confirm their suitability in pre-sterilized systems
- Materials satisfy USP Class VI criteria
- Fittings meet industry requirements
- Aseptic disconnect operation

SPECIFICATIONS


- **Materials:**
 - USP Class VI, ADCF Polysulfone
 - Platinum-cured, USP Class VI silicone o-rings
 - 316 stainless steel springs
- **Leak Resistance:**
 - 125 psi, 8.62 bar
- **Operating Temperature:**
 - -40° F (-40° C) to 280° F (138° C)
- **Sterilization:**
 - Gamma radiation up to 50kGy, autoclavable up to 270° F (132° C) for a maximum of 60 minutes

Non-Valved Connectors

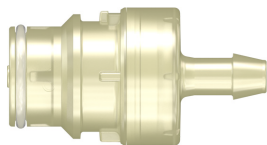


SDAM655NPS-001

Male Open Flow Connector, 600 Series Barb, 1/4" (6.4 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

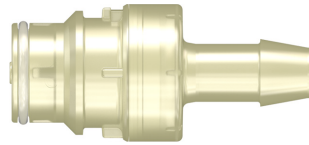

SDAM670NPS-001


Male Open Flow Connector, 600 Series Barb, 3/8" (9.5 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

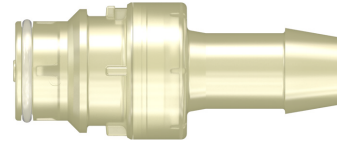

SDAM680NPS-001


Male Open Flow Connector, 600 Series Barb, 1/2" (12.7 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

Valved Connectors

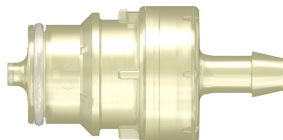


SDAVM655NPS-001

Male Valved Connector, 600 Series Barb, 1/4" (6.4 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

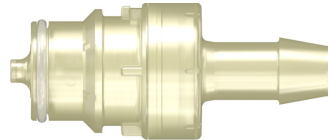


SDAVM670NPS-001

Male Valved Connector, 600 Series Barb, 3/8" (9.5 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

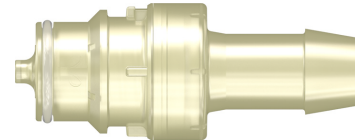

SDAVM680NPS-001

Male Valved Connector, 600 Series Barb, 1/2" (12.7 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

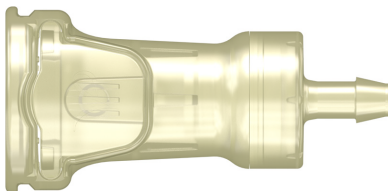


SDAVF655NPS-001

Female Valved Connector, 600 Series Barb, 1/4" (6.4 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAVF670NPS-001


Female Valved Connector, 600 Series Barb, 3/8" (9.5 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

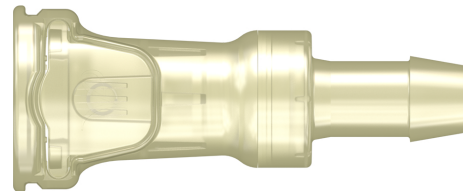

SDAVF680NPS-001


Female Valved Connector, 600 Series Barb, 1/2" (12.7 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring

Cap Connector

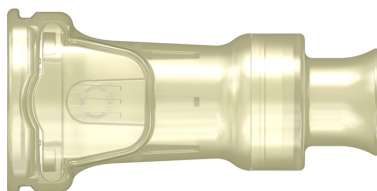

SDAFCNPS-001

Female Cap Connector, Animal-Free Polysulfone

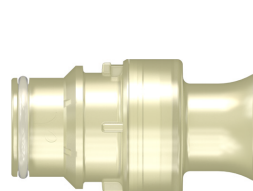

SDAMPNPS-001

Male Plug Connector, Animal-Free Polysulfone, Silicone O-Ring

Plug Connector